



WAX-310 WAX HEATER

Model: **WAX-310**

Certificate No.: **PTC24112705101E-EM01**

TECHNICAL DATA:


Power: **100W**

Rated input:

AC 110-240V 50/60Hz 1A

Capacity: **500 ml**

Temperature Range: **55-140°C**

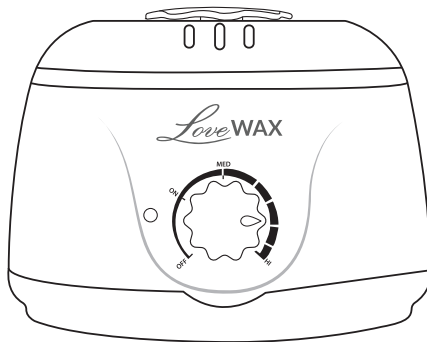

IMPORTANT SAFETY INFORMATION:

When using electrical appliances, especially in the presence of children, always follow basic safety precautions.

Read the entire instruction manual before using the device.

1. Do not place or store the appliance where it may fall or be pulled into a bathtub, sink, or any other container filled with water.
2. Do not immerse the appliance in water or any other liquid. Keep the appliance and surrounding area dry at all times.
3. If the appliance falls into water, unplug it immediately. Do not reach into the water to retrieve the appliance.
4. As with all electrical appliances, electrical components may remain live even when the appliance is switched off. Always use caution during operation.

WARNINGS:

To reduce the risk of burns, electric shock, fire, or injury:

1. Do not use the appliance if:
 - the power cord or plug is damaged,
 - the appliance is not operating properly,
 - the appliance has been dropped or damaged,
 - the appliance has come into contact with water.
2. Close supervision is required when the appliance is used near children.
3. Keep the power cord away from heated surfaces.
4. Do not immerse the device in water or other liquids.
5. Do not allow the power cord to hang over the edge of a table or countertop to prevent accidental pulling or tripping.
6. Operate the device only on a stable and level surface.
7. After use, disconnect the device from the power supply. Allow it to cool before cleaning or storing.
8. Do not wrap the power cord around the device during operation.
9. Do not use accessories not recommended by the manufacturer.
10. Always test the wax temperature before use.
11. Use the device only for its intended purpose as described in this instruction manual.

Operating Instructions:

1. Place the protective metal collar onto the rim of the wax heater.
2. Remove the plastic lid from the device.
3. Insert the wax container into the heater using the handle for easier handling.
4. Set the temperature control knob to the maximum setting until the wax reaches a liquid consistency. The melting process takes approximately 15 minutes. To accelerate heating, close the lid during operation.
5. Once the wax has completely melted, adjust the temperature control to a medium setting to maintain the proper consistency throughout use.
6. Before starting the treatment, test the wax temperature by applying a small amount to the inside of your wrist. Adjust the temperature if necessary.
7. After use, remove the wax container using the handle.
8. Switch the device off by turning the control knob to the OFF position. Replace the lid to protect the wax from dust and contaminants.
9. Clean the device only after disconnecting it from the power supply and allowing it to partially cool down.
10. Before the next use, repeat steps 1–7.

CAUTION:

Hot wax may cause burns.

- For professional use only.
- Read the instruction manual carefully before use.
- Do not leave the device unattended while in operation.



In accordance with art. 13 sec. 1 of the Act of 11 September 2015 on used electrical and electronic equipment (Journal of Laws 2015 item 1688), we inform you about the correct handling of waste electrical and electronic equipment:

1. It is prohibited to place used electrical and electronic equipment together with other waste - this is confirmed by the marking in the form of a "crossed-out bin", which requires selective collection of this type of waste.
2. Electrical and electronic devices may contain hazardous substances, mixtures and components, which, after entering the environment, may pose a serious threat to the health and life of humans and living organisms.
3. Used electrical and electronic equipment should be transferred only to authorized collection points, the list of which should be included on the website of each Municipal Office.
4. Households play a key role in the waste management system for used electrical and electronic equipment due to the possibility of their direct transfer to authorized collection points and the elimination of undesirable social habits resulting in leaving waste equipment in places not intended for this purpose.

IMPORTER:
GnBLAB sp. z o.o.
ul. Piotrkowska 270
90-361 Łódź, Polska
hurt.gnb-lab.com
uwagi@gnb-lab.com
Made in PRC